



More than just tubing, the disposable ambIT® Cassette is the pump mechanism of the ambIT® Infusion System

Volumetric accuracy of +/- 6% ensures safe and accurate infusion flow rates.



The small increments of medication that move through the tubing, and around the rollers in the cassette, allow for a very precise amount of medication to be dispensed.

Safety Features:

Disk mechanism prevents free flow

Cassette prongs ensure that there is only one way to install cassette and tubing

Allows pharmacies to use validated USP 797 processes, reducing infection risks and increasing patient safety when filling sterile medication solutions.

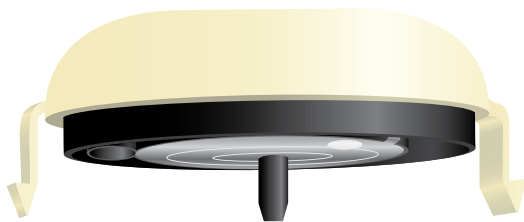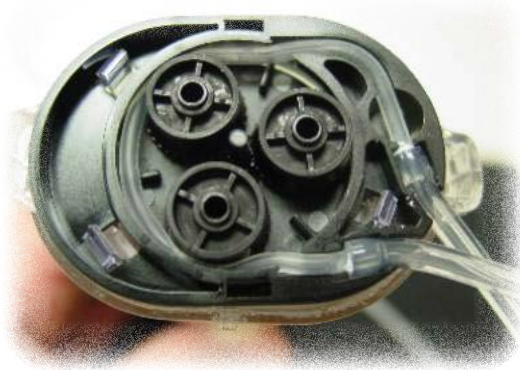

Flexibility and Convenience:

- U.S. and internationally patented technology for detecting occlusions or other pump malfunctions
- Accommodates standard and non-standard collapsible fluid medication bags
- Anti-free flow protection
- Quick manual priming
- Optional built-in air elimination filter
- Latex-free and DEHP-free fluid path

To be used in conjunction with our full product line of ambIT® electronic ambulatory infusion pumps.

504 West 8360 South
Sandy, Utah 84070 | USA
801-352-1888
877-352-1888
Fax: 801-352-1818
service@ambitpump.com

ambIT® is a division of
Summit Medical Products, Inc.



Available Cassette Configurations

ambIT® Cassette, Filterless, Spike
ambIT® Cassette, Filter, Spike
ambIT® Cassette, Filter, Spike with Yellow Striped Tubing
ambIT® Cassette, Filterless, Male Luer
ambIT® Cassette, Filter, Male Luer
ambIT® Cassette, Filter, Spike with Blue Striped Tubing
Luer to Spike Adapter

Specifications

| | |
|----------------------------|--|
| Cassette Tubing Length: | 60 inches (152 cm) to patient 10 inches (25 cm) to medication bag |
| Volume within Tubing line: | 2.3 ml |
| Cassette Tubing Material: | PVC with polyethylene liner |

ambIT® Cassettes are manufactured in accordance with internationally recognized volumetric accuracy standards (IEC 60601-2-24), unlike elastomeric and mechanical pumps.

This device is restricted to sale by or on the order of a physician.

504 West 8360 South
Sandy, Utah 84070 | USA
801-352-1888
877-352-1888

Fax: 801-352-1818
service@ambitpump.com

ambIT® is a division of
Summit Medical Products, Inc. 